









## **3 SECTION ELECTRIC EXAMINATION COUCH** MODEL: AM-03




This electric examination table has a very strong table designed to withstand constant heavy use with minimal maintenance. This examination table addresses the requirements of a wide range of handling issues to benefit both patient and doctor.

## **Technical Specifications:**

Product Size	196*66*50/89 cm		
Packing Size	199*74*53 cm		
CBM/pc	0.8		
Weight/kgs	N.W.:75 KGS G.W.:90 KGS		
Motor QTY	1 motors		
Adjustable Position	1 for up & down		
Bed surface	3 sections		
Power plug	Customize		
Voltage	220V		
Color of leather	Blue		
Material of leather	PU/PVC leather, Synthetic leather		



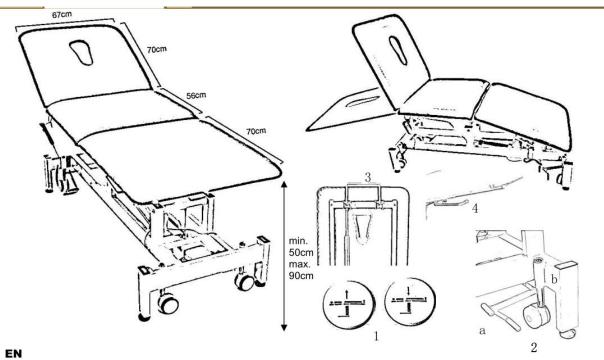

- \* The height up and down control by motor remote padel (1),the arrow up means up
  \* The backrest up and down control by hydraulic lever (3)
  \* The end section up and down control by hydraulic lever (4)

- \* The foot lever(2) control wheels and down, push the 4 foot levers to position (a) can move the bed to other places; when using the bed, please keep them in position (b)